



IBP SETS

INVASIVE BLOOD PRESSURE MONITORING SYSTEMS

Monitoring of arterial blood pressure is a mainstay of hemodynamic monitoring in acutely or critically ill patients. Sky Medical offers full product portfolio for invasive blood pressure monitoring addressing any invasive blood pressure monitoring procedure in ICU, OR and during interventional cardiology.

SKyMedical



Color coded tubing

Color coded stopcock

Ergonomic easy-to-use cable connection

- Integrated pressure transducer with an easy to use flush device and 3-way color coded stopcock
- Transparent transducers housing for good visualization
- Flush device continuous flush rate of 3ml/h or 30ml/h at 300mmHg provides simple and safe priming

SET FEATURES & BENEFITS

Our sets provide continuous and accurate measurement of physiological pressure and determination of other hemodynamic parameters. Our pressure transducers are fully integrated with flush device and 3-way stopcock, our sets vary from single to triple channel transducers, color coded tubing, color coded stopcocks, needle-free access for blood sampling and closed systems with blood reservoir. Customized sets according to the clinical procedure and the hospital demands are available.



Integral reservoir for blood aspiration

Color coded stopcock

Ergonomic mounting on the transducer plate

Self-flushing needle free sampling port


Closed blood sampling sets

Closed blood sampling sets combine safe and efficient blood sampling with integrated aspiration syringe and needle-free sampling port




Precise and safe real-time blood pressure measurement to manage critically ill patients - in one device.


ORDERING INFORMATION

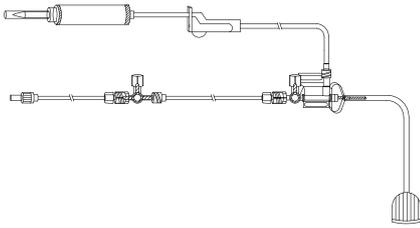


BP-S00003
Single transducer set

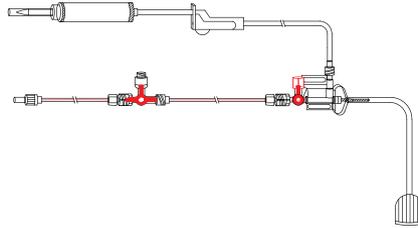


BP-S00028
Single transducer set with
self-flushable needle free valve

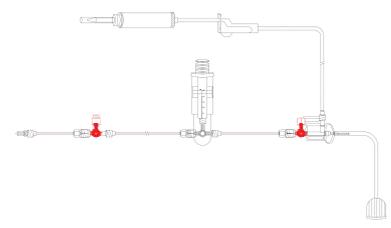

BP-S00031
Single transducer closed
blood sampling set

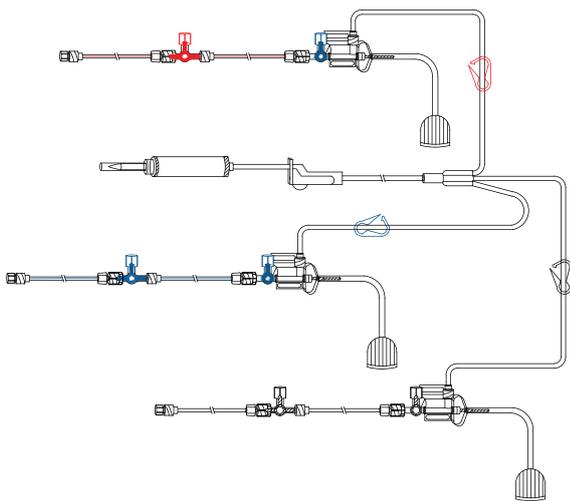


BP-T00001
Triple transducer set

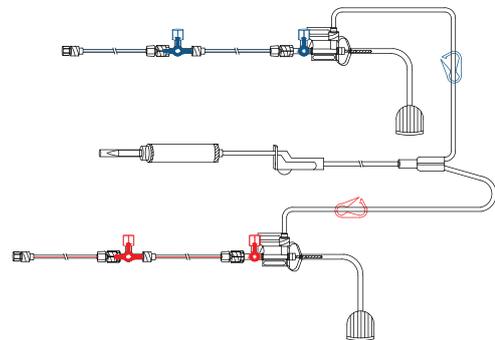


BP-D00001
Double transducer set



BP-S00032
Disposable pressure
transducer

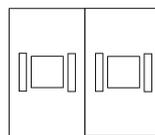

DO-D0002
2 plates transducer
organizer

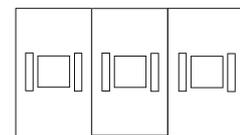

DO-T0003
3 plates transducer
organizer



We are offering more configurations to create efficient and safe procedure.

DO-H0000
Pole hanger for
transducer plate



CONTACT US:

Viceny 24, 931 01 Šamorín, Slovak Republic

Phone : +421 312 851 700 Email : info@med-sky.com Web: www.med-sky.com

SKyMedical
6BBP02 rev.2 07/2020